

COL

| REF | vials for kits | vial filling volume | drops per vial |
|--------|----------------|---------------------|----------------|
| COL10 | 10 | 0 - 1,45 ml | 45 |
| COL20 | 20 | 0 - 1,45 ml | 45 |
| COLB10 | 10 | 0 - 2,50 ml | 80 |
| COLB20 | 20 | 0 - 2,50 ml | 80 |

COL.C

| REF | vials for kits | vial filling volume | drops per vial |
|---------|----------------|---------------------|----------------|
| COLC30 | 30 | 0 - 1,45 ml | 45 |
| COLC40 | 40 | 0 - 1,45 ml | 45 |
| COLC50 | 50 | 0 - 1,45 ml | 45 |
| COLCB30 | 30 | 0 - 2,50 ml | 80 |
| COLCB40 | 40 | 0 - 2,50 ml | 80 |
| COLCB50 | 50 | 0 - 2,50 ml | 80 |



KIT FOR THE
PREPARATION
OF ALIQUOTS
OF EYE DROPS
WITH FRESH AUTOLOGOUS
HEMOCOMPONENTS

COL / COL.C

BIOMED DEVICE S.R.L.
via Ippolito Nievo, 89
41124 Modena - Italy

phone +39.059.343929 - fax +39.059.4550010
e.mail: info@biomeddevice.it
www.biomeddevice.it




PRACTICAL QUICK SAFE COL / COL.C



Resealable vials

up to
20
vials

standard drops
0,031
ml

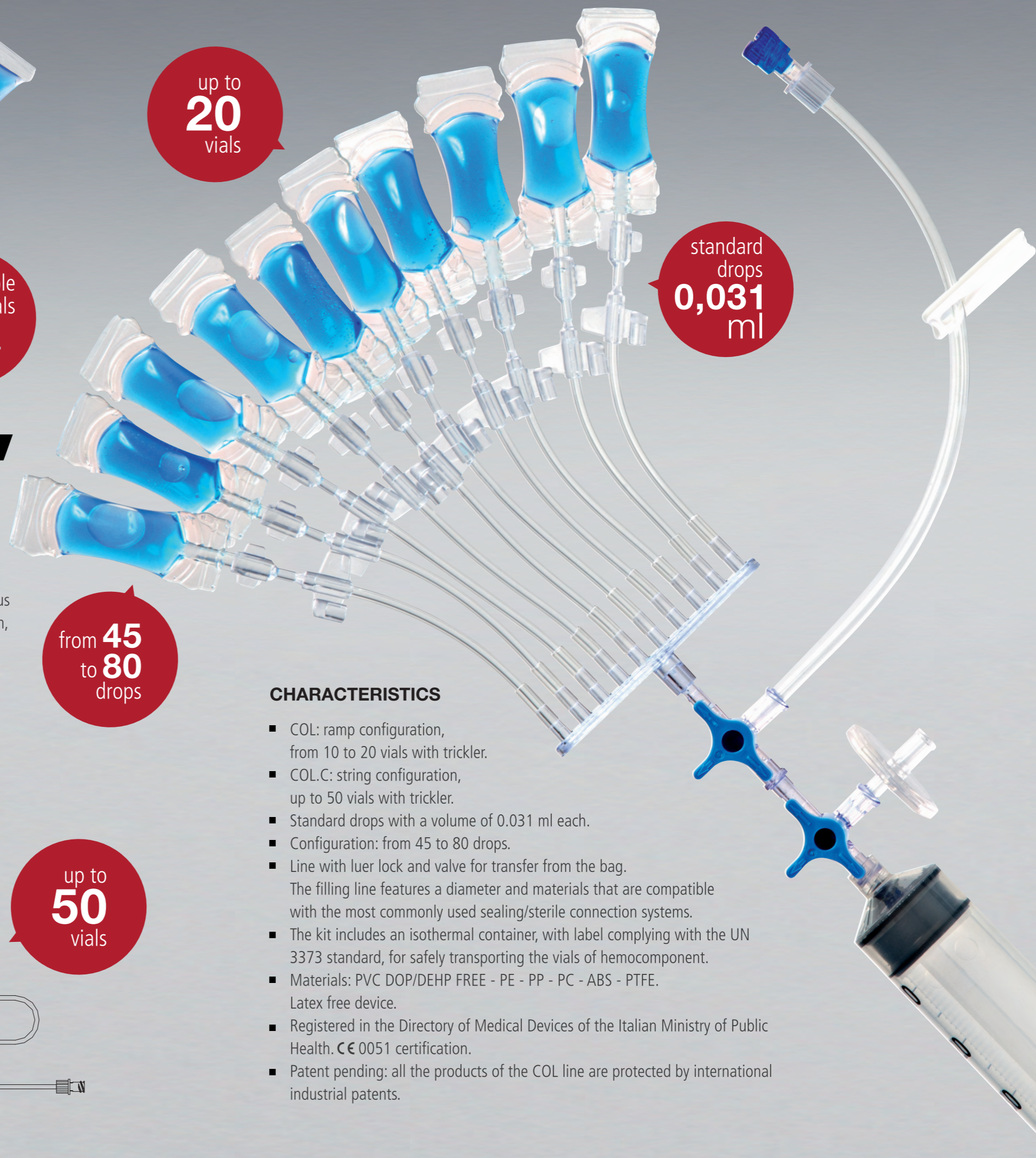

from **45**
to **80**
drops

up to
50
vials

KIT FOR THE PREPARATION OF ALIQUOTS OF EYE DROPS WITH FRESH AUTOLOGOUS HEMOCOMPONENTS. SYSTEM FOR HEMOCOMPONENTS PREPARATION AND PRESERVATION.

COL is the result of the experience and long-term research activities of Biomed Device. COL allows for the preparation of aliquots of eye drops using fresh autologous hemocomponents (serum, platelet-rich plasma (PRP), platelet concentrate, autoserum, serum eye drops). It is a closed-circuit, pre-connected and sterile kit that is very practical to use. **Each technical detail has been especially designed to fulfil the requirements of hospital workers in the most simple, fast and safe manner.**

CHARACTERISTICS

- COL: ramp configuration, from 10 to 20 vials with trickler.
- COL.C: string configuration, up to 50 vials with trickler.
- Standard drops with a volume of 0.031 ml each.
- Configuration: from 45 to 80 drops.
- Line with luer lock and valve for transfer from the bag. The filling line features a diameter and materials that are compatible with the most commonly used sealing/sterile connection systems.
- The kit includes an isothermal container, with label complying with the UN 3373 standard, for safely transporting the vials of hemocomponent.
- Materials: PVC DOP/DEHP FREE - PE - PP - PC - ABS - PTFE. Latex free device.
- Registered in the Directory of Medical Devices of the Italian Ministry of Public Health. CE 0051 certification.
- Patent pending: all the products of the COL line are protected by international industrial patents.

DO YOU NEED MORE VIALS?

Choose the COL.C configuration to obtain up to 50 vials!

